




CFDA CE

BEIEN

Focus on Minimally Invasive Surgery

ENDO-ARTISAN™

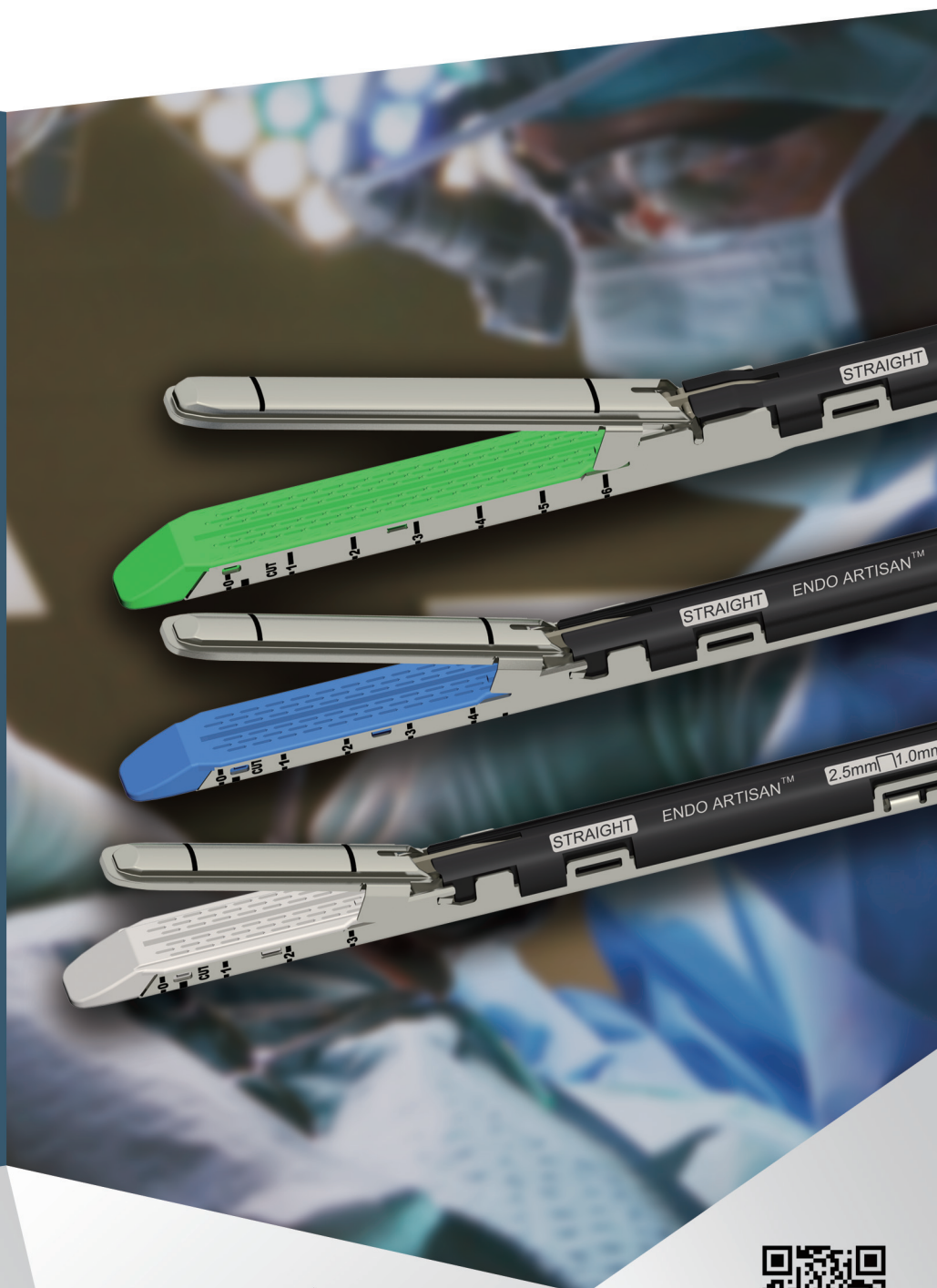

Disposable
Endoscopic Reload Unit
Straight

Contact us

+0086-510-8378-2500

+0086-510-8378-2511


Service@Beien-Surgery.com

Wuxi Beien Surgery Device Co.,Ltd.

*Quality First
Deliver life-affirming
Result*



www.beien-surgery.com

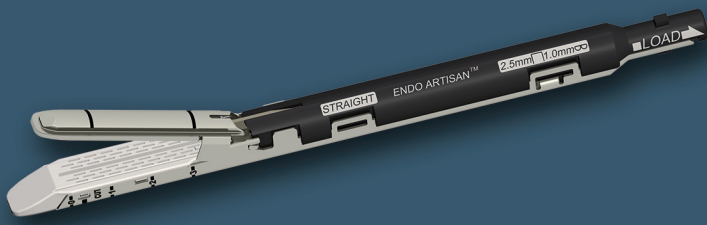

PRODUCT INFORMATION


RELOAD SELECTION CHART

THIN RANGE	MEDIUM RANGE	THICK RANGE
VASCULAR		
2.5 mm	3.5 mm	4.8 mm

12mm port

15mm port

- 2.5 : Artery, Vein, Venous plexus, Ligament, Omentum, Mesentery etc.
- 3.5 : Stomach, Small intestine, Colon, Anastomosis, Bronchus, Lung parenchyma etc.
- 4.8 : Distal rectum, Pyloricus, Trachea, Lung parenchyma etc.

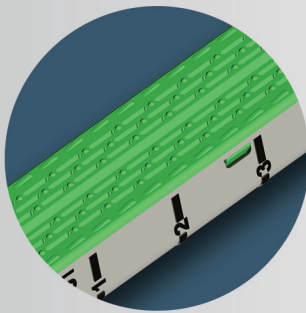

Convex Surface,
Effectively Anti-Slip

FEATURES

- Independent blade in each reload to provide the sharpness and ensure complete tissue cutting, reduce tissue damage.
- 6 rows of staple lines for both tissue and vascular sulus to maximize the safety of anastomosis and closure.
- Titanium alloy staplers, bigger tension-resistance on anastomosis.

SPECIFICATION

- Suture Length: 30-60mm
- Trocar Size: 12mm, 15mm
- Compressed Tissue Size: 1.0-2.0mm
- Staple Quantity: 6 rows; 48pcs, 66pcs, 90pcs.

Model	Code	Description	Cartridge Color	Stapler Size
WS30	3025L	30mm Straight Reload Unit	White	2.5mm
WS45	4525L	45mm Straight Reload Unit	White	2.5mm
WS60	6025L	60mm Straight Reload Unit	White	2.5mm
BS30	3035L	30mm Straight Reload Unit	Blue	3.5mm
BS45	4535L	45mm Straight Reload Unit	Blue	3.5mm
BS60	6035L	60mm Straight Reload Unit	Blue	3.5mm
GS45	4548L	45mm Straight Reload Unit	Green	4.8mm
GS60	6048L	60mm Straight Reload Unit	Green	4.8mm